




FORA G71

Blood Glucose Monitoring System

Owner's Booklet
ver 2.0


APPEARANCE AND KEY FUNCTIONS OF THE METER

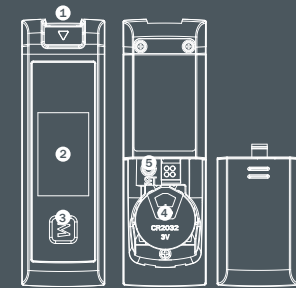

- 1 Test Slot**
Is where you insert the test strip. The meter will turn on automatically after insertion.
- 2 Display Screen**
Guides you through the test using symbols and simple messages.
- 3 M Button**
Located at the front of the monitor and marked with an "M", is used to turn on the monitor and enter the memory.
- 4 Battery**
Replace the battery when the power is low.
- 5 Set Button**
Is used to set up the meter.

NOTE
The meter will turn off automatically after 120 seconds without any action or you can press and hold the main button for 3 seconds to turn off the meter.

DISPLAY SCREEN

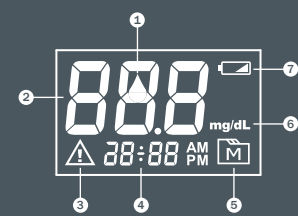


- 1 Blood Drop Symbol**
Flashes when meter is ready for sample application.
- 2 Test Result Area**
Test results are displayed here.
- 3 Error Warning**
Appears with error messages.
- 4 Month / Day / Hour / Minute**
- 5 Memory**
Indicates a test result stored in memory.
- 6 Measurement Unit**
The measurement unit is fixed as mg/dL.
- 7 Battery Symbol**
Appears when the battery is low.

FORA G71 TEST STRIPS

Your system measures the amount of sugar (glucose) in whole blood. Blood is applied to the absorbent hole on the test strip and is automatically drawn into the reaction cell where the reaction takes place.

The test strip consists of the following parts:

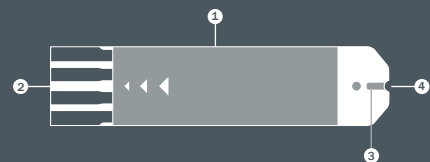

- 1 Test Strip Handle**
Hold this part to insert the test strip into the slot.
- 2 Contact Bars**
Insert this end of the test strip into the meter. Push it in firmly until it will go no further.
- 3 Confirmation Window**
This is where you confirm if enough blood has been applied to the absorbent hole in the strip.
- 4 Absorbent Hole**
Apply a drop of blood here. The blood will be automatically absorbed.

FORA • Blood Glucose Meter • Five-Year Warranty Registration Card

Complete this warranty card and mail it promptly

Name _____
Street Address _____
APT# _____ City _____ Country _____
State _____ Zip _____ Home Telephone _____
Date of Purchase (Month/Day/Year) _____
Meter Serial Number (See label on back of meter) _____

I would like to receive special offers and updates from FORA

www.foracare.ch

FORA Suisse AG

ForaCare Suisse AG

Neugasse 55, 9000 St. Gallen, Switzerland

PLACE
STAMP
HERE

TESTING YOUR BLOOD


Before testing, be sure to read this section and the test strip insert found in the test strip box carefully. Make sure you have all the items needed for the test.

Test Steps

Step 1 Insert the Test Strip

Insert the test strip into the test port with the contact bars end first and facing up. The bars must be fully inserted into the port for an accurate result. The meter turns on automatically and displays the blood drop symbol.

LCD full display → flashing "▲"

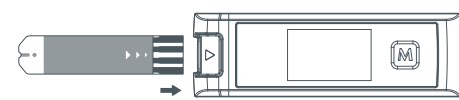

Step 2 Apply Blood Sample

Obtain a drop of blood of at least 0.5 µL using a lancing device. To obtain blood from other sites other than the fingertip, use the clear cap provided and refer to the strip package insert for more details.

When the "▲" symbol appears on the display, touch the drop of blood with the absorbent hole of the test strip. Do not rub the test strip with your finger or try to apply a smeared sample.

The blood should completely fill the confirmation window before the meter begins to count down. If the confirmation window is not filled completely before the meter begins to count down, do not add more blood to the test strip. Discard the test strip and start again. If you have trouble filling the test strip, please call your local customer service number for assistance.




NOTE

If you do not apply a blood sample within 3 minutes, the meter will turn itself off. You must remove the test strip and insert it back into the meter to restart the test.

Step 3 Get Results

The result of your blood glucose test will appear after the meter counts to 0. Blood glucose test results are automatically stored in the meter's memory.

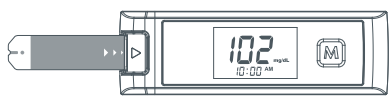

NOTE

Please do not change your treatment based on the result without first consulting your healthcare provider. Turn the meter off by removing the test strip. Discard the used test strip and lancet carefully according to your local regulations.

USING THE METER MEMORY


Viewing Results

Your Meter stores the latest 20 blood glucose test results with date and time in its memory. You can review the test results easily by the following steps.



Step 1 Enter the Memory Mode

With the meter turned off, press the main button twice. The first test result will appear, indicating that you are in the memory mode. If you continue to press the main button, you can then review the last 20 tests in the memory.

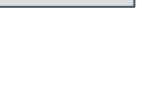


When using the meter for the first time or if the results have been deleted, "—" will appear, indicating that there are no test results in the memory.



Step 2 Recall Test Results

After the last test result, the most recent test result appears with the date and time. Press the main button to review your last 20 test results in order. When the memory is full, the oldest result will be deleted when the newest is added.

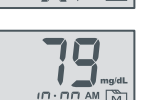

Step 3 Exit The Memory Mode

Press the main button for three seconds to turn off the meter.



NOTE

The results of control solution tests are NOT stored in the memory. Stored results are blood glucose results only.

IMPORTANT SAFETY INSTRUCTIONS

Read Before Use

Dear owner of a FORA G71 Blood Glucose Monitoring System,

The system consists of three main products: the meter, test strips, and control solutions. These products have been designed, tested, and proven to work together as a system to produce accurate blood glucose test results. Only use FORA G71 test strips and control solutions with the FORA G71 blood glucose monitoring system.

Intended Use

This system is intended for external use (in vitro diagnostic use) only. It is used for the quantitative measurement of glucose in samples of fresh capillary whole blood taken from the palm, forearm, upper arm, calf, thigh, or fingers. It is not intended to diagnose or screen for diabetes mellitus, or to be used on neonates. This system provides you with plasma equivalent results.

The Standard Kit for This System Includes

- 1 Meter
- 2 Owner's Booklet
- 3 Storage Case
- 4 Lancets
- 5 Warranty Card
- 6 Control Solutions (optional)
- 7 Test Strips
- 8 Lancing Device with a Clear Cap

NOTE

Check your system to ensure that it is unopened prior to use and that it contains all parts listed above. If anything is missing or damaged, please return your system to the place of purchase.

CHECKING THE SYSTEM WITH FORA CONTROL SOLUTIONS


FORA Control Solution contains a known amount of glucose that reacts with test strips and is used to ensure your meter and test strips are working together correctly.

Do a control solution test when:


- you first receive the meter,
- at least once a week to routinely check the meter and test strips,
- you begin using a new vial of test strips,
- you suspect the meter or test strips are not working properly,
- your blood glucose test results are not consistent with how you feel, or if you think the results are not accurate,
- practicing the testing process, or
- you have dropped or think you may have damaged the meter.

How to Perform a Control Solution Test

Step 1 Insert Test Strip
Insert a test strip into the test port with the contact bars end first and facing up. (Contact bars must be inserted all the way into the meter or you may get an inaccurate test result.) The meter will turn on automatically and display the followings in sequence:
LCD full display → flashing "▲"

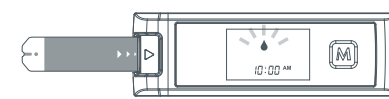


When the "▲" symbol appears on the display, press the main button and "CTL" will appear on the display. When the "CTL" sign is displayed, the meter will not store your test result in the memory. If you decide not to perform a control solution test, press the main button again and the "CTL" sign will disappear.

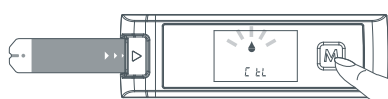

NOTE

Every time you perform a control solution test you must enter "CTL" test mode so that the test result will not be stored in the meter memory.

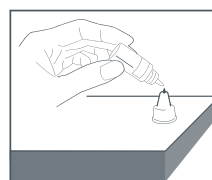

Step 2 Apply Control Solution

Shake the control solution vial well. Remove the cap from the control solution bottle. Place the cap on a flat surface. Squeeze the bottle, discard the first drop, and wipe the dispenser tip clean to prevent contamination. Squeeze out another drop and place it on the tip of the cap.

Hold the meter to touch the drop with the test strip's absorbent hole and the drop will automatically drawn into the test strip. Make sure the confirmation window is completely filled. The meter begins to count down.

NOTE

To avoid contaminating the control solution with the content of the test strip, it is recommended that you discard the first drop and place another drop of control solution on a clean surface first. Then touch the drop with the test strip.

Step 3 Check Results

After the meter counts down to 0, the control solution test result will appear. Compare the result with the range printed on the test strip vial. The result should fall within this range.

SPECIFICATIONS

Dimension: 86.9 (L) x 27.0 (W) x 11.3 (H) mm

Weight: 19 g (excluding battery)

Power Source: One CR2032 3V Lithium battery (for at least 500 measurements)

Display: LCD

Memory: 20 measurement results with date and time

Automatic detection of electrode insertion
Automatic reaction time count-down
Temperature warning

Operating Conditions: 10°C to 40°C (50°F to 104°F), below 85% R.H.
Storage/Transport Conditions: -20°C to 60°C (-4°F to 140°F), below 95% R.H.
Measurement Units: mg/dL
Measurement Range: 20-600 mg/dL (1.1-33.3 mmol/L)

When the battery is low, you can still do up to 50 tests.

The device has been tested to meet the requirements of ISO 15197, and electrical and safety requirements of: IEC/EN 61010-1, IEC/EN 61010-2-101, EN 61326-1, EN 61326-2-6.

ForaCare Suisse AG
Neugasse 55, 9000 St.Gallen, Switzerland
TEL: +41-71-220-1001
FAX: +41-71-220-1075
www.foracare.ch

CE 0459 IVD For self-testing 311-4274100-021

SYMBOL INFORMATION

SYMBOL	REFERENT
IVD	For in vitro diagnostic use only
Do not reuse	
Read instructions before use	
Keep away from sunlight	
Keep dry	
Temperature limitation	
Use by / Expiry date	
Use within 3 months after opening	
Batch code	
Manufacturer	
Serial number	
Caution, consult accompanying documents	
Sterilized using irradiation	
CE mark	
Do not use if package is damaged	
Green Dot: German dual waste disposal system	

SETTING THE METER AND DELETING THE MEMORY

Start with the meter off. Then press the setting button located on the back. The meter is now in the setting mode.



Step 1 Set the Year

The year will appear first, with the year flashing. Press and release the main button to add one year. To go faster, press the main button down longer. Once the correct year is displayed, press the setting button and the date will appear on the display with the month segment flashing.



Step 2 Set the Month

Press and release the main button until the correct month appears. To go faster, press the main button down longer. Once the correct month is displayed, press the button and the day will flash.




Step 3 Set the Date

Press and release the main button until the correct day appears. To go faster, press the main button down longer. Once the correct day is displayed, press the setting button and the time will appear on the display with the hour flashing. (To set the time, you must first set the date.)



Step 4 Select 12h or 24h Clock

Press and release the main button to select the 12h or 24h clock. These are different time presentation formats. If you select 12h, AM and PM will be shown alongside the time. If you select 24h, the time format will be from 00:00 to 23:59 without showing AM or PM.

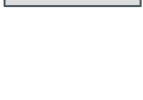

After you have finished selecting the 12h or 24h clock, press the setting button and then the hour will flash on the screen.



Step 5 Set the Hour

Press and release the main button to add one hour. To go faster, press the main button down longer. Once the correct hour is displayed, press the setting button and the minute will flash.



Step 6 Set the Minute

Press and release the main button to add one minute. To go faster, press the main button down longer. Once the correct minute is displayed, press the setting button to enter the memory deleting mode.



Step 7 Delete Results and Exit the Setting Mode

With the "no" symbol flashing on the display, you can easily delete all the results stored in the memory by pressing the main button once. The meter will display "yes". Then press the setting button to delete all the results and "no" will be displayed. The meter will then automatically turn off.

If you do not want to delete the results, do not press the main button but press the setting button again to turn off the meter.

When the meter is in the setting mode, it will also turn off automatically after being idle for two minutes.

DISPLAY MESSAGES AND PROBLEM-SOLVING GUIDE

The following is a summary of the display messages. If your meter displays an error message, please follow the actions for the error message as described in the table below. If the problem persists, please contact your local customer service agent for help.

MESSAGE	WHAT IT MEANS	WHAT TO DO
E-E Δ	Problem with the meter.	Review the instructions and try again with a new test strip. If the problem persists, please contact your local customer service agent for help.
E-D Δ	You may have removed the strip after applying blood to the absorbent hole.	Keep the test strip in the meter after it has absorbed the blood until the meter displays the test result.
E-F Δ	Memory error.	Please contact your local customer service agent for help.
E-U Δ	A used test strip is inserted, or there is an electrical problem.	Repeat test with a new test strip. If the error message re-appears, please contact your local customer service agent for help.
E-L Δ	The temperature of the environment, meter, or test strip is outside the system's operating range. You cannot perform a test until the meter and test strip are within the operating range of 10°C to 40°C (50°F to 104°F).	Repeat the test after the meter and test strip are within the operating temperatures.
E-L Δ	Result is below the measurement limit of 20 mg/dL (1.1 mmol/L).	This indicates hypoglycemia (low blood glucose.) Please seek medical assistance immediately.
H Δ	Result is above the measurement limit of 600 mg/dL (33.3 mmol/L).	This indicates severe hyperglycemia (high blood glucose.) Please seek medical assistance immediately.

For questionable results, please see test strip package insert for important information.