

# FORA<sup>®</sup> FORAFINE

Clinically  
Validated

## Pen Needles ACS022

### Enjoy Comfort!



### FORAFINE Pen Needles – Quality you can trust



- A universal design compatible with all major injection pen brands
- Thin Wall cannula for an optimized drug flow
- Color-coded for easy identification of various needle lengths
- A Double-layered Silicone Coating on the needle surface makes injecting extremely comfortable
- Available in 4 mm, 5 mm, 6 mm and 8 mm lengths





## Thin Wall Cannula

Improved medication through a greater flow rate

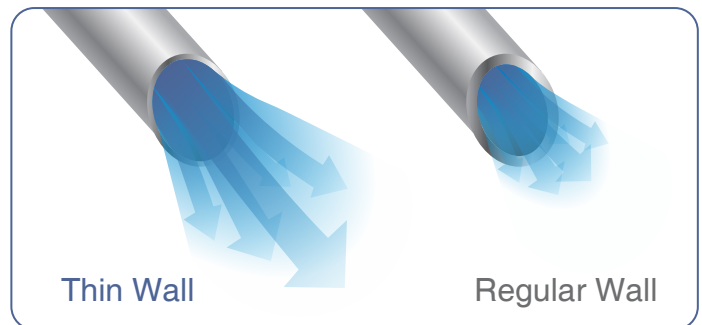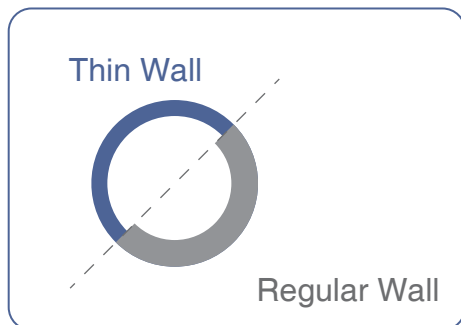

Advanced quality medical cannula:  
All 4 mm, 5mm, and 6 mm pen  
needles have a fine 32G cannula.


- Improved medication through a greater flow rate\*.
- Due to the large inner diameter of the cannula, less force is required when pushing the injection button.
- Supports the simple delivery of highly concentrated insulin doses of 200 Units/ml.

\*The flow rate test is performed according to Appendix A, ISO 11608-2.<sup>1</sup>




## Premium Quality and Features


Improve overall injection experience

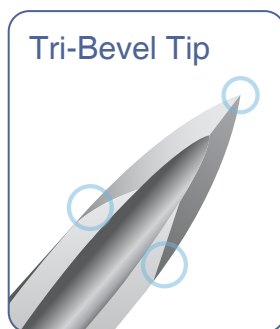

The Tri-bevel Tip and precise Metal Cutting of the needle greatly reduces the injection resistance.

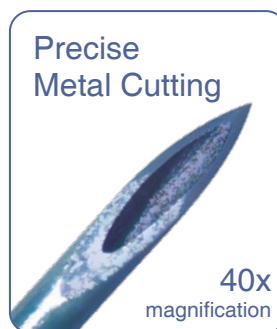

With its precise Metal Cutting and premium polish, the pen needle tip minimizes the injury to the skin.

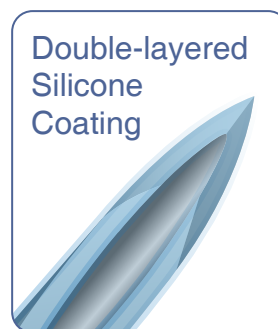

A Double-layered Silicone Coating on the needle surface offers an extremely comfortable injection experience.



The external protective cap with a tear-off tab allows the sealing peel-foil to be removed quickly and easily.

Reference:

1. International Organization for Standardization (2012). *ISO 11608-2:2012 - Needle-based injection systems for medical use — Requirements and test methods — Part 2: Needles*. Geneva, Switzerland: International Organization for Standardization.



## Safety and Reliable Quality

EN ISO13485:2016 and CE certified

The FORA® FORAFINE pen needles are EN ISO13485:2016 and CE certified.

The FORA® FORAFINE pen needles meet the ISO standards for compatibility. They are latex-free, non-toxic, and non-pyrogenic.



## Universal Twist

Fit all major brands of insulin pens

The FORA® FORAFINE pen needles are compatible with all major injection pen brands.

- Autopen®
- BD Pen®
- FlexPen®
- NovoPen® / NovoLet® / InnoLet®
- HumaPen Ergo® / HumaPen LUXURA® HD / HumaPen MEMOIR®
- Humalog KwikPen®
- SoloSTAR®
- SymlinPen®
- BerliPen®
- Omnican Pen®
- OptiClik® / OptiPen® / OptiSet®



## Choosing the correct pen needle length

Color-coded for easy identification

The FORA® FORAFINE pen needles are available in four different sizes: 4 mm, 5 mm, 6 mm, and 8 mm. The different pen needle lengths with their respective outer and inner diameters are color-coded for easy identification of the needle size.





## Product Overview

### Types and Product Range



| Items                     | Q'ty/Box   | Cannula Diameter   | Needle Tip |
|---------------------------|------------|--------------------|------------|
| FORA® FORAFINE 32G x 4 mm | 100 pieces | 0.23 mm, Thin Wall | Tri-bevel  |
| FORA® FORAFINE 32G x 6 mm | 100 pieces | 0.23 mm, Thin Wall | Tri-bevel  |
| FORA® FORAFINE 31G x 5 mm | 100 pieces | 0.25 mm, Thin Wall | Tri-bevel  |
| FORA® FORAFINE 31G x 8 mm | 100 pieces | 0.25 mm, Thin Wall | Tri-bevel  |

**ForaCare offers comprehensive Diabetes Care Solutions and services for people with diabetes. Foracare Full Care Solutions:  
Test Simply, Live better.**



Pen Needles




Diabetes Care



Fever Care



Hypertension Care

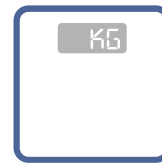

Weight Care



Telehealth Solution

## ForaCare Suisse AG

Neugasse 55, 9000 St. Gallen, Switzerland

TEL : +41-71-220-1001

FAX : +41-71-220-1075

URL : [www.foracare.ch](http://www.foracare.ch)

E-mail : [chservice@foracare.com](mailto:chservice@foracare.com)



[www.foracare.ch](http://www.foracare.ch)